


Brush that loosens cells in pancreatic cysts could become weapon against cancer

January 28 2021

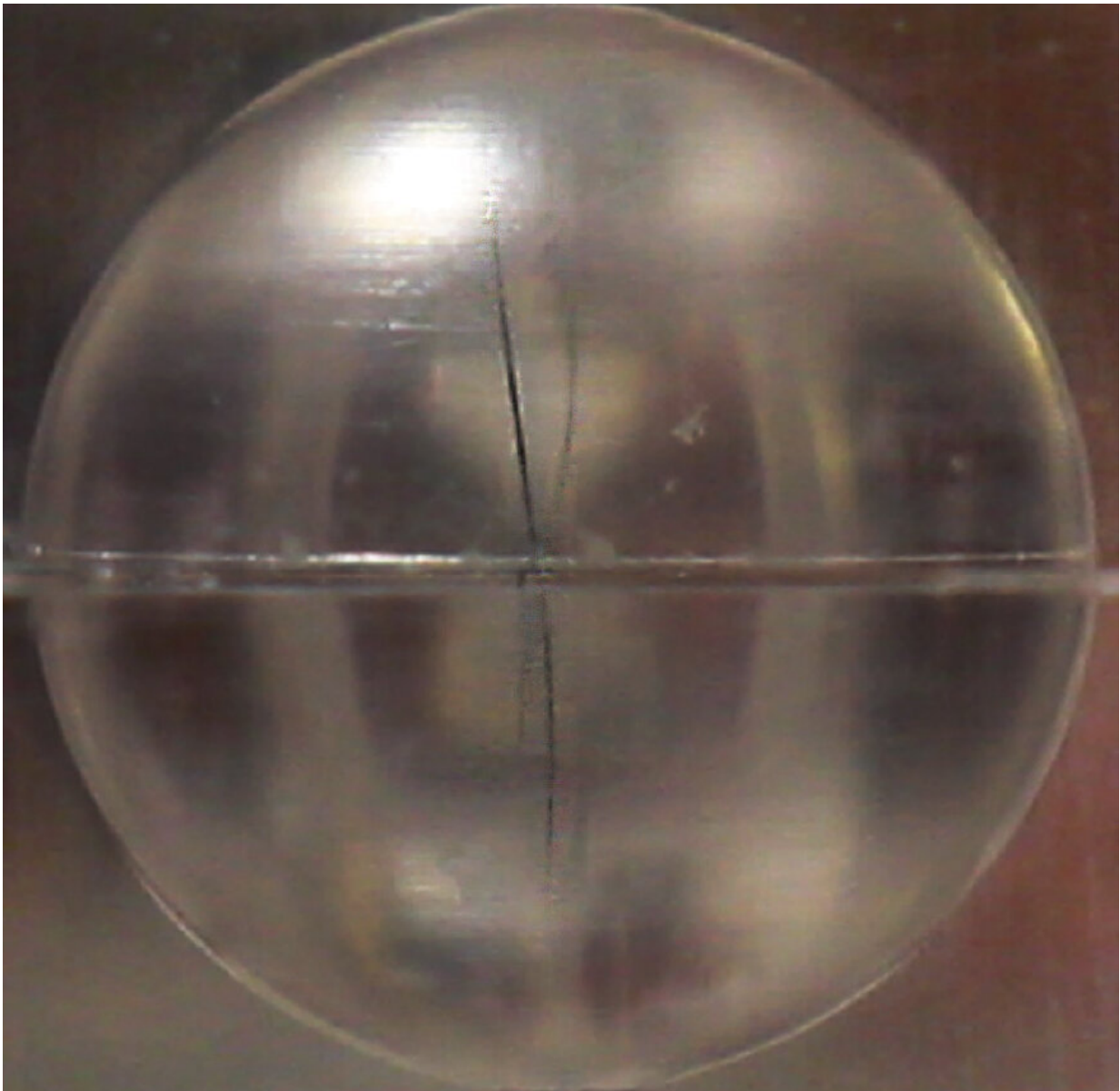

A model of a cyst that was used to test the loop brush. Credit: Filipe Marques

In a potential step forward for early pancreatic cancer diagnosis, a miniaturized brush was successfully tested as a possible technique for loosening cells from the walls of pancreatic cysts. If proven safe to use, the technique could dramatically increase the number of cystic cells available for analysis, a key advance given the scarcity of cells in cystic fluid drawn with a needle alone.

Researchers at Sweden's KTH Royal Institute of Technology proposed modifying endoscopicultrasound-guided [fine needle aspiration](#) with the addition of a loop brush which can be inserted into a cyst through the needle, says Wouter van der Wijngaart, professor of Micro and Nanosystems at KTH.


"In order to detect malignancy in a pancreatic cyst, the cells from the inner wall of the cyst need to be analyzed," van der Wijngaart says. "You may find some of those released into the cystic [fluid](#) naturally; but sometimes you find none—that is why we try to brush them off."

Pancreatic cancer is the fourth leading cause of cancer death in Western countries. With a survival rate of between 2 and 9 percent, early detection is critical.

"Pancreatic cystic lesions are precursors of pancreatic cancers and a prime target for early diagnosis," van der Wijngaart says. However, he says, 20 percent of cystic fluid samples contain negligible amounts of cells, or no cells at all, hampering a proper diagnosis.

To evaluate the functionality of the device, the researchers performed tests on models of cysts and intestines from pigs and cows. They fitted

the aspiration needle with a 50 micrometer-thin nitinol wire shaped in a 1 cm diameter loop, which is fixated to a 280 micrometer-thick nitinol guidewire. The loop is compressed and inserted into a cyst through a fine needle with a diameter less than 0.5 mm.

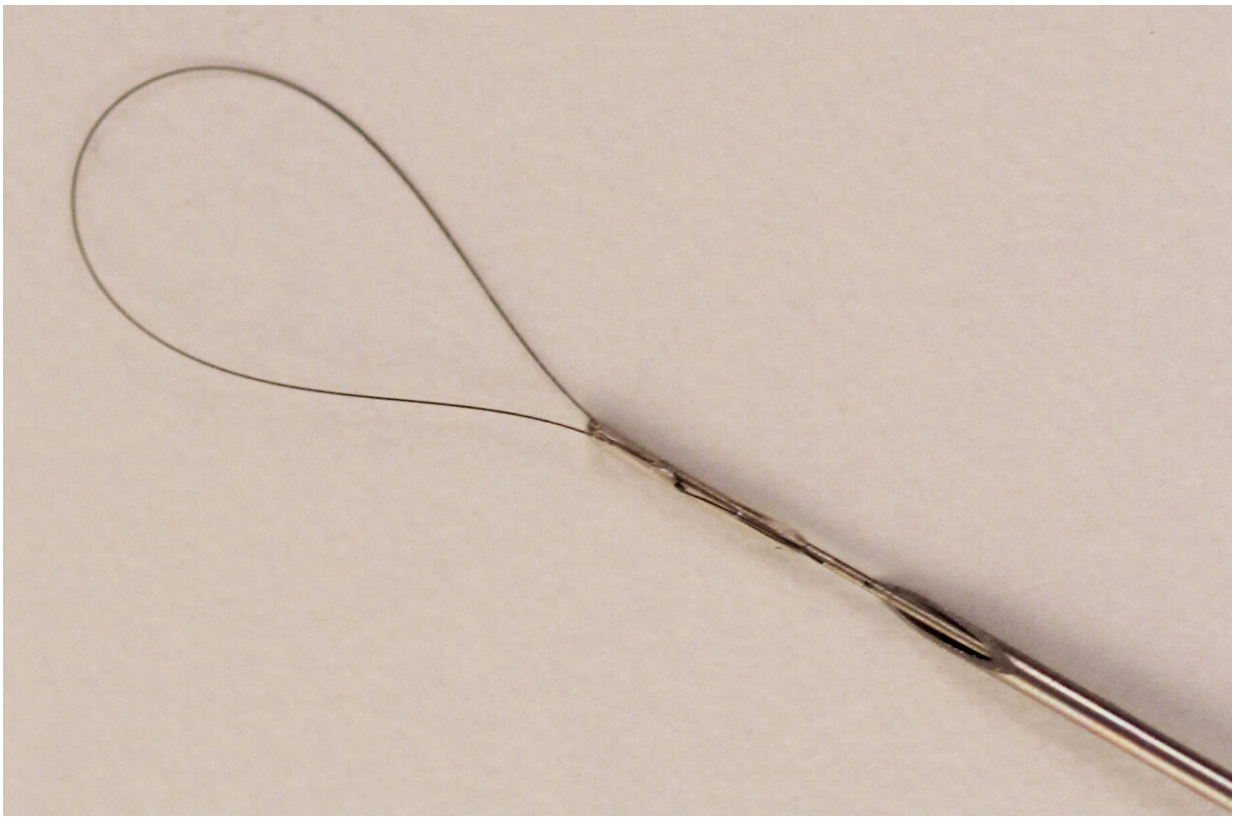

The loop brush extending from needle. Credit: Filipe Marques

Once inside the cyst, the brush opens up and is rotated with low force to abrade the [cells](#) from the cyst wall. Then the loop is retracted and the cystic fluid is aspirated through the [needle](#) and collected for cytology.

The cell content increased up to 54 times when using low-viscous water as fluid in soft tissue, soft-like cavities; up to 174 times when using

high-viscous glycerol; and at least 10 times in cow ovary cysts.

The project was initiated by co-author Francisco Baldaque Silva, surgeon at Karolinska University Hospital, and the experimental work was performed by Ph.D. student Filipe Marques. Also contributing to the study were KTH researcher Niclas Roxhed and Urban Arnelo, a researcher at Karolinska Institutet. The project was sponsored by Region Stockholm.

More information: Filipe Marques et al. A loop-shaped minimally invasive brush for improved cytology sampling of pancreatic cysts during EUS-FNA, *MEDICAL DEVICES & SENSORS* (2020). [DOI: 10.1002/mds3.10165](https://doi.org/10.1002/mds3.10165)

Provided by KTH Royal Institute of Technology

Citation: Brush that loosens cells in pancreatic cysts could become weapon against cancer (2021, January 28) retrieved 27 April 2024 from <https://medicalxpress.com/news/2021-01-loosens-cells-pancreatic-cysts-weapon.html>

<p>This document is subject to copyright. Apart from any fair dealing for the purpose of private study or research, no part may be reproduced without the written permission. The content is provided for information purposes only.</p>
--