


# A new microfluidic system could keep tabs on cancer treatment

August 3 2022, by Catherine Graham

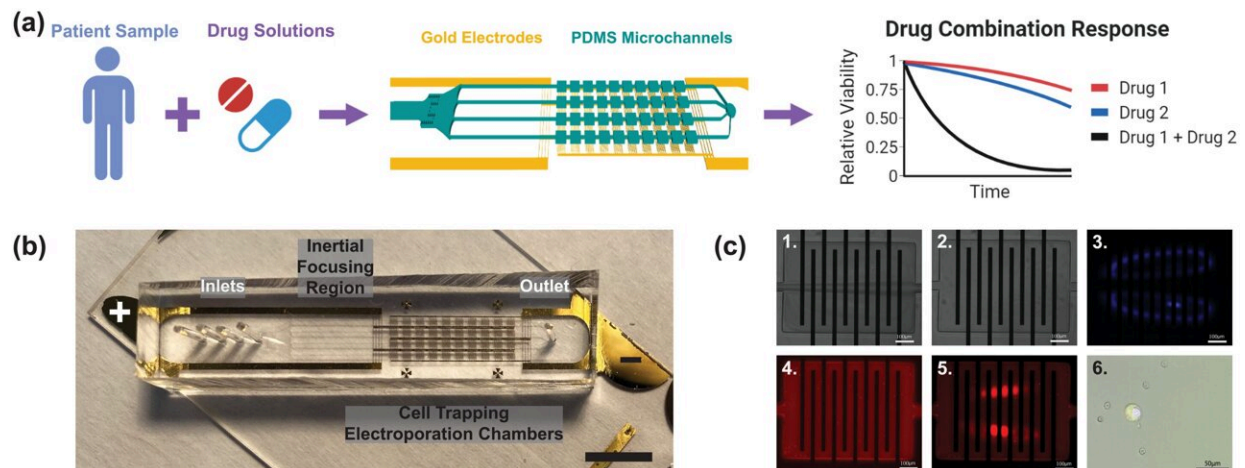

Overview of workflow for assessing drug response from cancer cells. (a) Schematic workflow for directly assaying diseased cells from patient liquid biopsy samples for subsequent drug combination screening. Cell-trapping PDMS microfluidic chambers (turquoise) and Au electrode arrays (gold) of the assembled device allow cell purification from blood and subsequent drug response analysis (b) Overhead photograph of the device. The device consists of an inlet region accommodating multiple solution ports with coarse filter for blocking debris, inertial focusing region for size-based focusing of larger cancer cells from smaller blood-borne contaminant cells, and cell trapping region (10x4 chamber arrangements) fitted with patterned microelectrodes (10 planar microarrays per chamber) for cytosolic delivery of bioreagents via electroporation. Scale bar represents 5 mm. (c) Representative time-lapse panel images of the integrated process consisting of 1. cell trapping, 2. flushing away blood cells, 3. visualization of vortices, 4. molecular incubation, 5. visual confirmation of successful electroporation, and 6. merged microscopic image of

both background blood cell (white arrow) and a viable, electroporated cancer cell (black arrow), indicated by overlapping red (Propidium Iodide) and green (Calcein AM) fluorescent signals. Scale bar represents 100  $\mu\text{m}$  for panels 1. - 5. and 50  $\mu\text{m}$  for panel 6. Credit: *PLOS ONE* (2022). DOI: 10.1371/journal.pone.0264907#sec020

Found in a patient's blood, circulating tumor cells, or CTCs, have powerful potential to indicate how well a cancer treatment is working and guide treatment decisions. Even so, CTC technologies have not yet been widely implemented in clinics, in part because CTCs are very rare, and separating them from the billion normal cells in blood is a challenge.

A new microfluidic system developed by Johns Hopkins engineers could make it easier to collect and test [drug](#)-resistant CTCs from blood samples known as "liquid biopsies," which are faster, less invasive, and safer than traditional tissue biopsies.

The team's results appeared in *PLOS ONE*.

"Our workflow has the potential to directly evaluate circulating [tumor cells](#) and determine new, optimum drug combinations that make [cancer cells](#) responsive to treatment again," said study principal investigator Soojung Claire Hur, the Clare Boothe Luce Assistant Professor in the Department of Mechanical Engineering at the Whiting School of Engineering.

Because they can be taken more frequently, liquid biopsies help clinicians evaluate patients at different phases of their [cancer treatment](#). For example, these biopsies can glean information on why some cancer [cells](#) become resistant to a patient's current treatment, said Hur.

Her lab had previously developed a technique to isolate CTCs for testing purposes: a [blood sample](#) flows through a microfluidic chip, which recognizes and traps only the rare tumor cells, like a panhandler sifting for gold.

In this new study, the researchers set out to adapt that technique to investigate cells' real-time drug response. First, the team ran whole blood samples spiked with a drug-resistant lung cancer cell line through the previously mentioned purification process. Then they used their new microfluidic device as an electroporator, creating holes in the cancer cell membranes with jolts of electricity before directly injecting the cells with a combination of two cancer drugs.

They found that combining the right drugs could overcome resistance and restore therapeutic response in the tumor cells, said Hyun Woo Sung, a Ph.D. candidate in Hur's lab and lead author of the study.

According to Sung, the system's key innovation is its reliance on a single-cell imaging pipeline, a technique that identifies discrete cell-to-cell differences within the cell population. In this study, they used this approach to determine whether every cell they collect and electroporate is permeable for drug delivery or not, and whether the cell is alive or dead after the jolt.

"Generally, a higher voltage will lead to more individual 'instance' of cells we see in our images that are permeabilized, but likewise a higher 'instance' of cells that die from too big of a jolt. By collecting all this information from the cells at a range of input voltages, our approach allows us to calculate the electroporation performance and identify the optimal voltage that work best for the cells in our system," Sung said.

The single-cell analysis is also promising in that it addresses the rare nature of CTCs in [blood samples](#). Conventional techniques for analysis

require a certain number of cells to be present in order to successfully process the sample. In contrast, the Hur lab's system offers an alternative way to analyze a smaller population of cells, such as CTCs.

Although this study provides a proof-of-concept, with more work the researchers believe the system will someday give clinicians the information they need to determine whether individual patients will respond to a specific drug or not.

"In future studies, we plan to increase the volume of blood processed per run and expand on the types of therapies we are delivering to the cells, beyond drugs. Our technology could also be used to monitor response to other treatments, like [gene therapy](#) and immunotherapy, just from a simple [blood](#) test," Sung said.

**More information:** Hyun Woo Sung et al, Sensitizing drug-resistant cancer cells from blood using microfluidic electroporator, *PLOS ONE* (2022). [DOI: 10.1371/journal.pone.0264907#sec020](https://doi.org/10.1371/journal.pone.0264907#sec020)

Provided by Johns Hopkins University

Citation: A new microfluidic system could keep tabs on cancer treatment (2022, August 3) retrieved 25 April 2024 from <https://medicalxpress.com/news/2022-08-microfluidic-tabs-cancer-treatment.html>

|  |
|--|
| <p>This document is subject to copyright. Apart from any fair dealing for the purpose of private study or research, no part may be reproduced without the written permission. The content is provided for information purposes only.</p> |
|--|